DIRECTIONS FOR USE

ALL-IN-ONE VACUUM MATTRESS

English: Directions for use. **Ferno All-in-One Vacuum Mattress**
Deutsch: Bedienungsanleitung. **Ferno All-in-One Vakuummatratze**
Français: Mode d'emploi. **Matelas à dépression Ferno All-in-One**
Español: Instrucciones de uso. **Colchón de vacío todo en uno Ferno**
Svenska: Bruksanvisning. **Ferno All-in-one vakuummadrass**
Italiano: Istruzioni per l'uso. **Ferno All-in-One Materassino A Depressione**
Dansk: Brugsanvisning. **Ferno, All-in-One vakuummadrass**
Nederlands: Gebruiksaanwijzing. **Ferno Alles-in-één Vacuüm Matras**
Suomi: Käyttöohje. **Ferno All-in-One -tyhjiöpöytä**
Norsk: Brukerveiledning. **Ferno All-in-One vakuum-madrass**

Distributor:

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English Directions for use

1. Introduction

The Ferno All-in-One Vacuum Mattress is a combined stretcher and vacuum mattress, hence the name All-in-One. The mattress can be on the stretcher at all times working as a stretcher mattress.

If vacuum immobilization is needed it is available immediately. No extra storage space for a vacuum mattress is required. The immobilization is well suited for fixation of patients when bone and spinal injuries are suspected. Deep tissue damage can also be avoided.

The vacuum mattress is fitted with 6 handles on both longer sides. It is also equipped with four belts with Velcro. Quick release buckles and/or handles at each short end are also available on request.

Two persons are able to carry the mattress a shorter distance and at longer distance additional carrying belts can be used.

You can pull the mattress on the ground by the handles or by connecting a carrying belt or rope.

The All-in-One Vacuum Mattress consists of an outer shell of airtight material that is filled with an inner construction containing granules. When evacuating the air out of the mattress, it can be shaped around the patient.

When the air has been evacuated from the All-in-One Vacuum Mattress, it is rigid and it is not possible to mould it any more. The mattress wraps around the patient offering an optimal support and stabilization.

The mattress can carry a patient in water. (See precautions section).

2. Precautions

- Never mould the All-in-One Vacuum Mattress on the top of the head or under the foot soles as it might cause undesirable pressure against the spine.
- On a current basis check that the vacuum pump works correctly and that it is possible to get the air out of the All-in-One Vacuum Mattress to make it hard.
- Avoid handling the mattress next to sharp objects. If the mattress is pulled be aware of obstacles on the ground, such as stones, sticks, glass, sharp edges etc.
- Check pump, All-in-One Vacuum Mattress and carrying straps as to damages and wear and tear. Repair damages. Change defect parts immediately.
- Avoid handling the mattress next to open fire.
- On a current basis check that the mattress has no cracks.
- The mattress is unstable in water. You can not leave the patient unattended in water.
- At transportation in Aircraft, please note that the atmospheric pressure can be lower and the mattress is not enough rigid. In this situation evacuate more air from the mattress with your pump.
- Avoid carrying a patient by using the handles at the short ends only. An undesirable bending might occur.

3. Description

The vacuum mattress consists of an outer shell made of hard-wearing PVC coated fabric and a unique granulate filled innerbag.

When the air is evacuated the granulate is stuffed into a tight substance making the mattress rigid and unable to shape. Before and while the air is evacuated the mattress may be moulded according to the shape of the patient. An optimum support of the patient is obtained when all air has been evacuated.

The mattress is equipped with 6 handles on both sides. It is also equipped with four color coded belts with Velcro. When not in use, park the belts and fasten them by use of the velcro in the gap between the row of handles and the mattress.

In case you do not want to use the stretcher mattress facility, Ferno has made a little unique feature that makes it easy to fold and pack the mattresses in a protective bag together with pump and repair kit so it only takes up little space.

4. Use of the All-in-One Vacuum Mattress

1. Put the mattress on an even surface near the patient.
Check that the valve is open.
2. Mould the mattress so that it is even or if necessary in an appropriate state for the situation.
3. Close the valve by turning it to the right (clockwise) and connect the pump.
4. Place the injured patient in the medically desired position on the mattress.
5. Mould the mattress around the patient so that it becomes a stable and comfortable unit to the patient.
Mould up some material from the mattress between the legs of the patient.
6. Close the Velcro straps according to the color code and tighten them. Chest belts are red and light blue.
Cross them over the chest.

Caution! Make sure that the material does not press on top of the head or under the foot soles.

7. Evacuate the air, and continue moulding the mattress. Tighten the belts once more.
Check that there is no additional pressure to the head or feet.
Check that the unit is stable and comfortable for the patient.
8. The patient is ready to be transported.

4.1 Use of the All-in-One Vacuum Mattress as a Stretcher/Transfer Mattress:

The All-in-One Vacuum Mattress can, when not evacuated, be used as a Stretcher/Transfer Mattress.
If it is the wide model:

1. Fold the longitudinal sides of the mattress, over the mattress so they meet in the middle.
2. Preferably place a sheet over the mattress.
3. Place the patient on the mattress and close the fastening straps from the stretcher.

If it is a narrower model simply place a sheet over the mattress.

5. Cleaning and Disinfecting

After each use the mattress should be wiped with a moist cloth and dried, if dirty the mattress may be washed with soapy water and rinsed with clean water.

Allow the mattress to dry completely before storage.

If the mattress needs disinfecting first clean the mattress carefully and let it dry.

Wipe the complete mattress with a solution containing 70% alcohol or use the disinfectant called LYSETOL AF.

6. Maintenance and Repair

At frequent recurrent intervals the vacuum mattress and the pump must be checked for visible damages.

At the same time the vacuum mattress must be tested as to tightness and the suction force of the suction pump must be checked according to Section 7.


Leaks can be found on the mattress by visual inspection or by reduced filling of the mattress with air, covering it with soapy water all over the surface and the welds watching out for air bubbles from the mattress. If filling of the mattress with air, this must be done without creating pressure on edges etc.

If leaks have been found on the mattress clean the area around the leakage and allow it to dry.

Find a patch in a suitable size and put glue on the patch as well as on the leakage on the mattress. Wait for the glue to become dust dry (approx. 2 minutes) and place the patch on the leakage. Firmly press the whole area of the patch establishing good contact and press all air out under the patch.

The repair will receive full strength after 24 hours.

Note: The belts are detachable, washable  and replaceable! If blood on the belts - put to soak in cold water first!

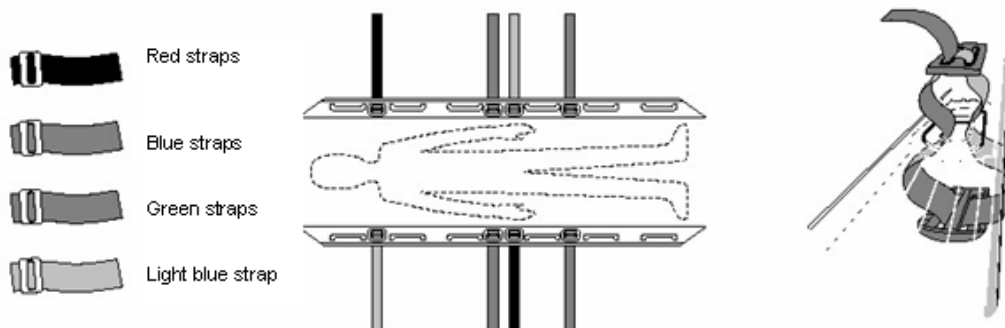7. Test of the All-in-One Vacuum Mattress and Pump

Straighten out the mattress, close the valve by turning it clockwise and connect the pump.

Evacuate the air until the mattress is rigid, leave it for about 2 hours. If the mattress is as rigid as it was when you started the mattress may be considered tight.

Open the valve by turning it counterclockwise.

If the mattress is not airtight see Section 6 Maintenance and Repair or check the valve for leakage.

Try to insert a new valve and repeat the test.

If the pump is not emptying the mattress check the pump by closing the connection with your fingertip and pump a few times. A negative pressure should be created in the tubing and can be felt at your fingertip. If the pump does not function do not repair it but return it for service to your dealer or to an Ferno service center.

8. Storage

After each use, check that the mattress is ready for next use. Loosen the valve by turning it counterclockwise and put the mattress on an even surface. When the granules have been distributed evenly, fold the mattress in its initial position.

Repack the mattress in order that it fits into your storage area.

The vacuum mattress is often positioned on the stretcher, or folded in a bag.

9. Specifications

The Ferno All-in-One Vacuum Mattress is tested in accordance with EN 1865

The Ferno All-in-One Vacuum Mattress is in conformity with Council Directive 93/42/EEC concerning Medical Devices

Dimensions:

Length x width x height: see product information below

Weight: see product information below

Heat resistance: 70 °C

Cold resistance: -30 °C

Capacity for foot pump: The pump has capacity to empty the mattress according to EN 1865

Materials:

Part	Material
Top and bottom material:	Polyester fibers PVC coated (flame resistant)
Extra reinforcement material:	Polyester fibers PVC coated (flame resistant)
Innerbag:	Polyamide

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