




AB GERMA®

Mini Transfer Stretcher/Mattress

We have together with our EMT customers developed a very handy, flexible and useful patient handling tool e.g. to move a patient from a stretcher to a bed or lift a patient from the floor bathtub or chair.

One big feature you will experience using Mini Transfer is when patient has to be rescued from a narrow environment e.g. bathroom. The great benefit is that you are able to bend the flexible stretcher/mattress around the patient and in step two lift and move the patient carrying in 10 options possible strong handles.

To make usage easy as possible, back of the Mini Transfer mattress is made with friction reducing material to slide easy behind or under the patient. Glass fibre material on the back is very durable and can be used sliding on floor or ground.

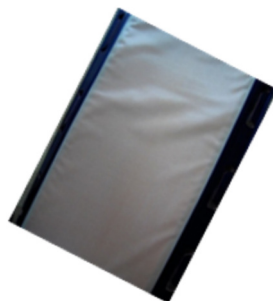

Glass fibre on the back.



Specification - Mini Transfer Stretcher/Mattress

Clear coloured green PVC-covered polyester (flame retardant)
Lie area is with 10 mm thick Polyether foam
Max loading capacity: 250 kg

Dimensions and weight

Length: 136 cm
Width: 78 cm
Weight: 1,8 kg

Ordering Specification

Mini Transfer Stretcher/Mattress

Art. No

271040000