



**AB GERMA®**

VACUUM MATTRESS  
AIO BLUE LINE



# PRODUCT SHEET

Germa Blue Line Vacuum Mattresses:

The Blue Line is Germa's Low-price range of Vacuum Mattresses. It is based on our AIO concept but with a less intricate inner bag, which requires a little adjustments of the granules before use. The lower price and less complex inner bag does still not tamper on Germa's high quality and long-life expectancy.

The intended use is as a Vacuum Mattress but to save storage space the mattress can lie on the stretcher at all times, acting as a stretcher mattress. If vacuum immobilisation is needed, the vacuum mattress is integrated in the stretcher mattress and is ready for use immediately.

## PRODUCT SPECIFICATIONS

Approved according to European Standard EN 1865 and EN 1789  
MDR 2017/745 Class 1 Medical device

### PRODUCT INFO

- 6 handles on each side – carry capacity 250 Kg
- Improved head and neck part for stronger immobilisation
- Improved chambers system

Warranty period: 3 years

### MATERIAL

- Outer material: Polyester fibre, PVC-coated, flame retardant
- Inner bag: Polyester
- Granules: Polystyrene granules (flame retardant)


### SPECIFICATIONS

- Length: 210 cm
- Width: 115 cm
- Height: 5,5 cm
- Weight 11 kg with restraint system
- Loading capacity: 250 kg

### TEST SPECIFICATIONS

- Maximum shrinkage on length: 1%
- Maximum shrinkage on width: 3%
- Working temperature of the mattress is -36°C to +75°C

# VACUUM MATTRESS AIO BLUE LINE

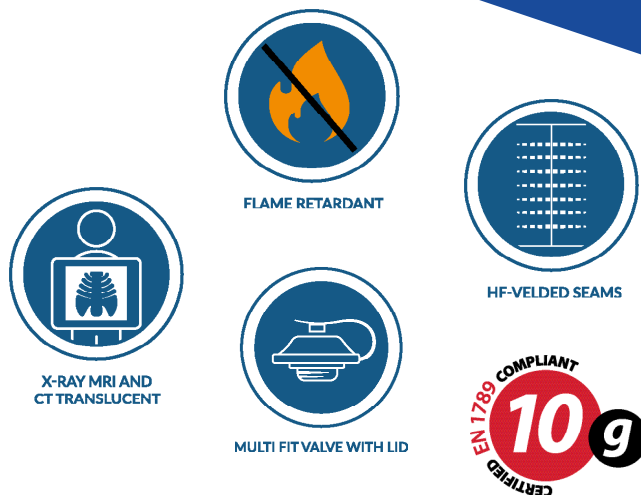

### ORDER INFORMATION

Art no: 271702000  
Vacuum Mattress AIO 100 Cm Blue Line

Art no: 271101001  
Vacuum Mattress AIO 100 Cm Blue Line with Valve in foot end

Art no: 271702001  
Vacuum Mattress AIO 100 CM Blue Line R2G Kit – including restraint system, footrest Vacuum Pump, repair kit, bag  
Weight: 13,5 kg

### SPARE PARTS

Art no: 12633800000  
Art: Germa vacuum valve Black

Art no: 271000501  
Art: Repair Kit Mattresses/Splints

Visiting address  
AB Germa  
Industrig. 54-56  
SE-291 36 Kristianstad, Sweden

Post address  
AB Germa  
Box 2127  
SE-291 02 Kristianstad, Sweden

www.germa.se  
Phone: +46 (0)44 12 30 30  
info@germa.se

## BENEFITS

Made to correspondent to EN 1789 and 1865  
Approved for European standard and crash test for 10-G

One-piece back, providing the patient with even spine support. Also allows preliminary X-ray diagnostics, direct in to x-ray no extra handling of patient.

Gives rescuers quick handling and excellent immobilisation for the patient.

Deliver consistent and rigid support to the patient  
Safe and convenient immobilisation.

Handles feature an ergonomic grip, making patient lifts as comfortable as possible for rescuers.  
Handles are reinforced so a patient can be lifted in a patient lifting device.

Handles fixed to the bottom of mattress, maintain immobilisation when lifting patient.

Robust welding—fully welded without any open seams or textiles, making hygienic cleaning of the mattress easy.

Valve with special designed hole, ensuring almost all suction devices fit directly into the valve, no adapters are needed. Reverse air lock system, no closing of valve at stressful site. Made from plastic makes it ideal for X-ray.

Safety belt fixation close to the mattress, safe and comfortable for patient.

PVC coated fabric makes it easy to clean and disinfect with most common chemicals.

Outer bag is made of flame retardant PVC-coated fabric. It is flexible making it easy to mould the mattress around the patient.

## VACUUM MATTRESS AIO BLUE LINE

Robust and efficient material to hold vacuum for more than 24 hours, last longer and safer for patient make trustfully immobilisation.

High load capacity 250 kg, safe to handle bariatric patients.

Two different closing mechanisms on the straps, either Velcro or plastic buckle, allowing the user to choose according to preference.

Repair Kit easy to use for effective repair, no return needed.

Folds easily, in case the stretcher mattress facility is not required.

High quality footrest aluminium vacuum pump evacuate with more than – 700 hpa, quick and time saving.