



AB GERMA®

VACUUM MATTRESS EASYFIX CHILD



PRODUCT SHEET

Germa offer new design Child EasyFix Vacuum mattress. Child EasyFix is developed from our earlier experiences from the original high quality materials traditionally used for producing Germa world class vacuum mattresses. Child EasyFix vacuum mattress using body shape design to optimize the vacuum mattress capabilities and limited storage volume. Intended immobilisation use for a child within the length of 90 cm up to 160cm this translate to a child between the age of approx. 5 years to 14 years. Child EasyFix unique construction which deliver efficient fixation of upper back, neck and head. Reinforced with base plate that cover cervical spine part. Together with fixation restraint system using X-system to support and create most effective immobilisation.


VACUUM MATTRESS EASYFIX CHILD



REINFORCED WOOD
INNER SYSTEM



FLAME RETARDANT

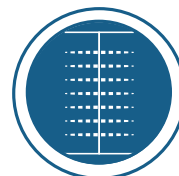

HF-VELDED SEAMS



X-RAY MRI AND
CT TRANSLUCENT



MULTI FIT VALVE WITH LID



PRODUCT SPECIFICATIONS

Approved according to European Standard
EN 1865 and EN 1789
MDR 2017/745 Class I Medical device

PRODUCT INFO

- 3 handles on each side – carry capacity 150 Kg
- Head strap system
- Colour coded Safety belt system with Fix lock – plastic buckles
- Head, neck and upper back part improved with wood part to secure stronger immobilisation
- Improved chambers system
- Light weighted and small storage

Warranty period: 3 years

MATERIAL

- Outer material: Polyester fibre, PVC-coated (flame retardant)
- Inner bag: Polyester
- Granules: Styropor granules (flame retardant)
- Wood board upper part: Birch plywood

SPECIFICATIONS

- Loading capacity: 150 kg
- Length: 160 cm
- Width: 85/48 cm
- Height: 4,5 cm
- 4,7 kg with restraint system

TEST SPECIFICATIONS

- Maximum shrinkage on length: 1%
- Maximum shrinkage on width: 3%
- Working temperature of the mattress is -36°C to +75°C

ORDER INFORMATION

Art no: 271504001

Vacuum Mattress Child EasyFix Ready2go

Weight: 7 kg

Art no: 271504000

Child EasyFix Vacuum Mattress

– including restraint system

SPARE PARTS

Art no: 12 633 800 000

Art: Germa vacuum valve Black

Art no: 15 000 000 000

Art: Easyfix repair kit

ACCESSORIES

Art no: 15 000 000 045

Art: Evacuation Straps, black, adjustment in length

Visiting address
AB Germa
Industrig. 54-56
SE-291 36 Kristianstad, Sweden

Post address
AB Germa
Box 2127
SE-291 02 Kristianstad, Sweden

www.germa.se
Phone: +46 (0)44 12 30 30
info@germa.se

BENEFITS

Made to correspondent to EN 1789 and 1865
Approved for European standard and crash test for 10-G

Germa head pillow design including rigid wood part solution.

Efficient, quick and improved neck and head immobilisation, with visible ears in one go.

Reinforced head, neck and upper back for fixation/immobilization of a patient.

Support upper back, neck and head for safe immobilization with built in birch plywood.

Easy to use, providing the patient with even spine support through the continuum of care.

Quick to prepare, allows preliminary X-ray diagnostics, direct in to x-ray minimize handling.

Germa improved chamber technology for quick and easy preparation.

Give rescuers quick handling and excellent immobilisation for the patient ensuring the granulate is evenly dispersed to give the patient stable and optimal immobilisation.

Deliver consistent and rigid support and immobilisation of the patient over time with limited risk.

Safe, premium immobilisation and convenience of patient fixation.

Handles with ergonomic grip.

Making patient lifts as comfortable as possible for rescuers.

All handles are reinforced so a patient can be lifted in a patient lifting device. Handles fixed to the bottom of mattress on piece.

Maintain rigid immobilisation when lifting the patient.

Robust welding—fully welded without any open seams or textiles. Making hygienic cleaning of the mattress easy.

Valve with special designed hole, ensuring almost all suction devices fit directly into the valve, no adapters are needed. Reverse air lock system, no closing of valve at stressful site.

Made from plastic makes it ideal for X-ray.

X-designed restraint belt fixation close to the patient, support immobilisation process and comfortable for patient. Supplied with four colour coded straps either with fix lock.

Helping the rescuer to place the straps correctly and easy to stretch and adjust.

VACUUM MATTRESS EASYFIX CHILD

Outer material is made of flame retardant PVC-coated fabric. It is safe and flexible making it easy to mould the mattress around the patient.

Design of the upper part allow the rescuer to place the arm(s) in two different options.

Make sure to give the rescuer access to patient arm for treatment requested via patient arm(s).

Robust and efficient material to hold vacuum for more than 24 hours.

Safety for patient and make confiding immobilisation.

High load capacity 150 kg, safe to handle patient.

Closing mechanisms of the straps, plastic buckle Fix lock. Allowing the user to quick tighten and adjust restraints on the patient.

PVC coated fabric makes it Easy to clean and disinfect with most common disinfectants

Repair Kit easy to use for effective repair, quick to repair no return needed.

High quality footrest aluminium vacuum pump evacuate with more than -700 hpa, durable and quick and time saving.

Visiting address
AB Germa
Industrig. 54-56
SE-291 36 Kristianstad, Sweden

Post address
AB Germa
Box 2127
SE-291 02 Kristianstad, Sweden

www.germa.se
Phone: +46 (0)44 12 30 30
info@germa.se



AB GERMA®

CHILD EASYFIX VACUUM MATTRESS

