




AB GERMA®

VACUUM MATTRESS EASYFIX LARGE

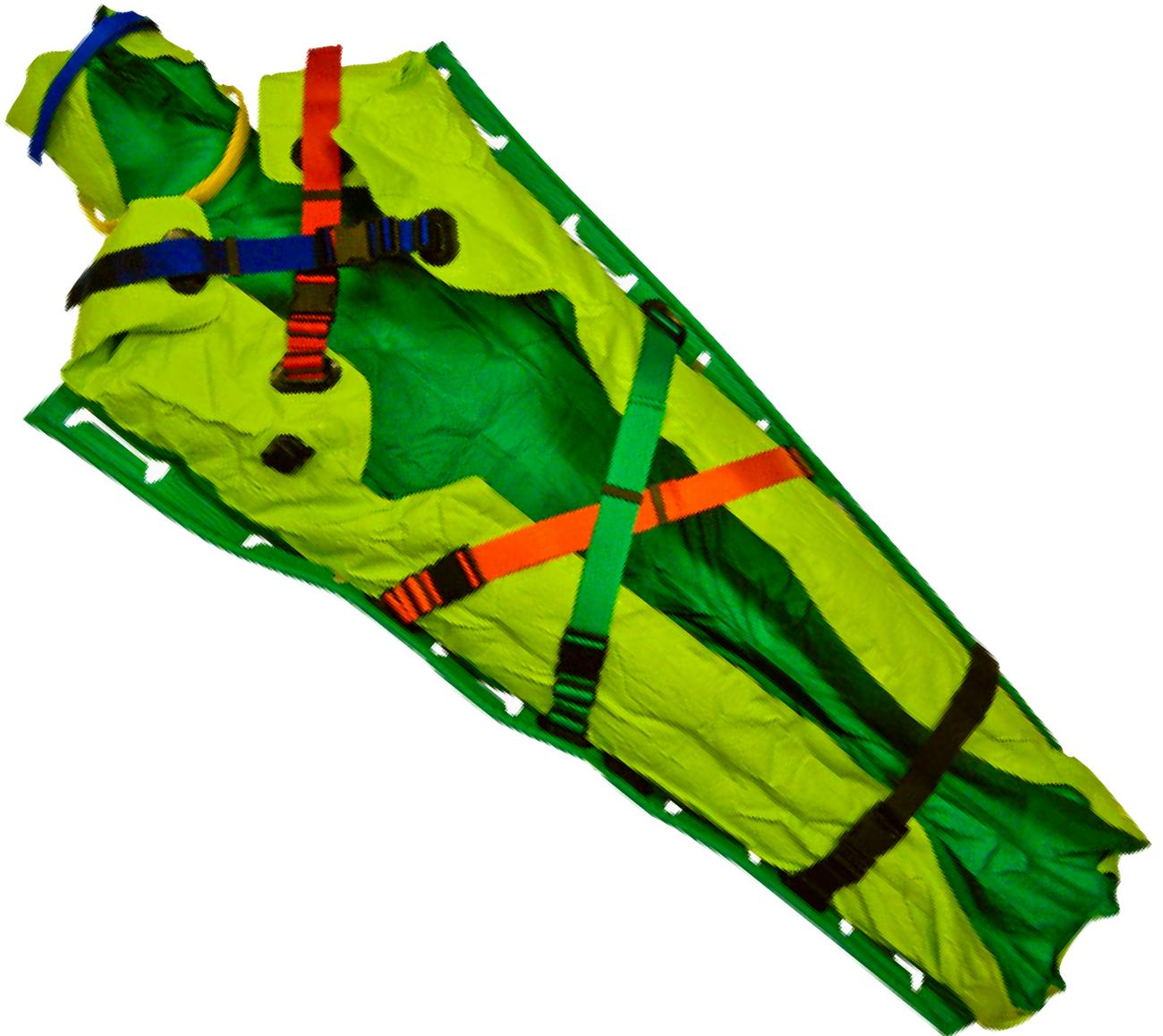

PRODUCT SHEET

Germa innovative renewed designed line of Vacuum Mattresses:

EasyFix is a development based on earlier experiences from the original high quality materials traditionally used producing Germa world class vacuum mattresses. EasyFix vacuum mattress using new body shape design to optimize the vacuum mattress capabilities. Offers 20% saving on storage space of the mattress compared to Germa All In One Vacuum Mattresses.

EasyFix unique construction which deliver improved fixation of upper back, neck and head by reinforced base plate that cover all of cervical spine. This together with fixation restraint system using x-system to support and create most effective immobilization.

PRODUCT SPECIFICATIONS

Approved according to European Standard EN 1865 and EN 1789
MDR 2017/745 Class 1 Medical device

PRODUCT INFO

- Warranty period: 3 years
- Head, neck and upper back part improved with wood part to secure stronger immobilisation
- New improved chambers system
- Compliant to EN 1789 and 1865
- Carry capacity of 250 kg, 10 handles
- Light weighted and small storage

MATERIAL

- Outer material:
- Polyester fibre, PVC-coated (flame retardant)
- Inner bag: Polyester
- Granules: Styropor granules (flame retardant)
- Wood board upper part: Birch plywood

SPECIFICATIONS SIZE:

- Length: 208 cm
- Width: 130 / 72 cm
- Thickness: 4,5 cm

VACUUM MATTRESS EASYFIX LARGE



REINFORCED WOOD
INNER SYSTEM



FLAME RETARDANT



HF-WELDED SEAMS



X-RAY MRI AND
CT TRANSLUCENT



MULTI FIT VALVE WITH LID



WEIGHT

- 8,0 kg without fixation straps
- 8,5 kg with fixation straps
- 9,6 kg total weight including fixation straps + footrest pump

TEST SPECIFICATIONS

- Maximum shrinkage on length: 1%
- Maximum shrinkage on width: 3%
- Working temperature of the mattress is -36°C to +75°C

ORDER INFORMATION:

Art no: 271501000

Art: Vacuum Mattress EasyFix Large

Art no: 271501001

Art: Vacuum Mattress EasyFix Large Ready2Go Kit
Including: Vacuum footrest pump+Repair Kit+ Bag

SPARE PARTS

Art no: 12 633 800 000

Art: Germa vacuum valve Black

Art no: 15000000000

Art: Repair Kit for EasyFix Vacuum Mattress

Visiting address
AB Germa
Industrig. 54-56
SE-291 36 Kristianstad, Sweden

Postal address
AB Germa
Box 2127
SE-291 02 Kristianstad, Sweden

www.germa.se
Phone: +46 (0)44 12 30 30
info@germa.se

BENEFITS

VACUUM MATTRESS EASYFIX LARGE

Approved according to European Standard
EN 1865 and EN 1789
MDR 2017/745 Class 1 Medical device

Germa head pillow design including rigid wood part solution.

Efficient, quick and improved neck and head immobilisation, with visible ears in one go.

One-piece back, providing the patient with even spine support.

Also allows preliminary X-ray diagnostics, direct in to x-ray no extra handling of patient.

Germa new chamber technology for quick and easy preparation. Gives rescuers quick handling and excellent immobilisation for the patient ensuring the granulate is evenly dispersed to give the patient stable and optimal immobilisation.

Ergonomic handles, making patient lifts as comfortable as possible for rescuers

Handles fixed to the bottom of mattress on piece, maintain rigid immobilisation when lifting the patient.

Robust welding—fully welded without any open seams or textiles, making hygienic cleaning of the mattress easy.

Valve with special designed hole, ensuring almost all suction devices fit directly into the valve, no adapters are needed. Reverse air lock system, no closing of valve at stressful site. As it is made from plastic makes it ideal for X-ray.

Improved V-shaped design to make sure upper wide part will make all sized patients the very best immobilisation of the back, neck and head.

X-designed restraint belt fixation close to the patient, support immobilisation process and comfortable for patient. Supplied with five colour coded straps with fix lock, helping the rescuer to place the straps correctly and makes it easy to stretch and adjust.

Outer material is made of flame retardant PVC-coated fabric. It is flexible making it easy to mould the mattress around the patient.

Design of the upper part allow the rescuer to place the arm(s) in two different options, making sure to give the rescuer access to patients arm for treatment requested via patient arm(s)

Extra layer built under the mattress, safety, minimize risk of damage the fixation, last longer.

Robust and efficient material to hold vacuum for more than 24 hours, safety for patient and make confining immobilisation.

High load capacity 250 kg, safe to handle heavy duty patients.

PVC coated fabric makes it easy to clean and disinfect with most common disinfectants, easy to keep clean and ready for use.

Folds easily, 20% less storage space required compared to All In One Vacuum Mattress, it is easy to fold and store the mattress in a protective bag and save storage space.

Visiting address
AB Germa
Industrig. 54-56
SE-291 36 Kristianstad, Sweden

Postal address
AB Germa
Box 2127
SE-291 02 Kristianstad, Sweden

www.germa.se
Phone: +46 (0)44 12 30 30
info@germa.se



AB GERMA®

VACUUM MATTRESS
EASYFIX LARGE

