




AB GERMA®

EASY TRANSFER MATTRESS

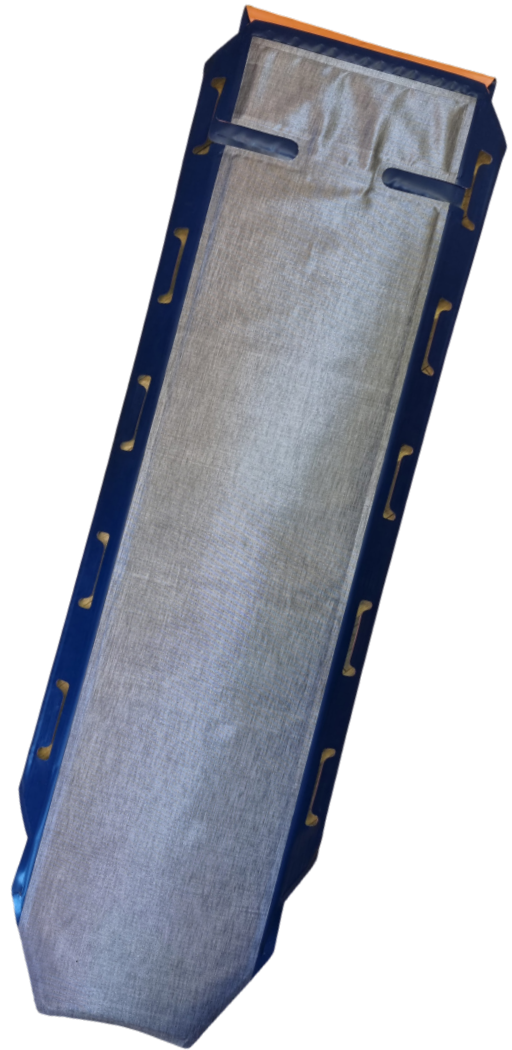

PRODUCT SHEET

Easy Transfer Mattress is designed in co-operation with professional users focusing on ergonomics patient handling, patient safety and comfort for patient.

Easy Transfer Mattress has been tested by trauma teams to secure efficient handling plus quick, safe and easy transfer of patients including heavy duty patients.

EASY TRANSFER MATTRESS



X-RAY MRI AND
CT TRANSLUCENT



FLAME RETARDANT



HF-VELDED SEAMS

PRODUCT SPECIFICATIONS

Approved according to European Standard EN 1865 and EN 1789
MDR 2017/745 Class 1 Medical device

PRODUCT INFO

- 5 handles on each side
- Glass fiber back side
- Carry capacity 250 Kg

Warranty period: 3 years

MATERIAL

- Outer material: Polyester fibre, PVC-coated, flame retardant.
- PP-plate and grass fibre material to reduce friction
- Inner material: 40 mm V - shaped Polyurethane Foam + 15 mm Sensus Foam

SPECIFICATION FOR THE UNI EASY TRANSFER MATTRESS

- Length: 185 cm
- Width: 50 cm
- Height: 6 cm
- Weight 4,7 kg

TEST SPECIFICATIONS

- Working temperature of the mattress is -36°C to +75°C

ORDER INFORMATION

Art no: 271030000
Easy Transfer Mattress

Visiting address
AB Germa
Industrig. 54-56
SE-291 36 Kristianstad, Sweden

Post address
AB Germa
Box 2127
SE-291 02 Kristianstad, Sweden

www.germa.se
Phone: +46 (0)44 12 30 30
info@germa.se

BENEFITS

EASY TRANSFER MATTRESS

Made to correspondent to EN 1789 and 1865
MDR 2017/745 Class 1 Medical devise

The PVC fabric is made strong and durable for long
and hard use.

The Uni Easy Transfer mattress has been
designed to fit most common stretchers
available on the market.

High level of patient comfort due to the
V-shaped foam plus the Sensus memory
foam.

5 handles on each side, allow paramedics
to easy lift, carry and slide the mattress with
patient on it.

For easier patient access during transport the
handle part fold under the mattress on the
stretcher.

Sliding surface on the reverse of fiberglass
material, reducing friction and allow sliding of
the mattress e.g. on the floor/ground or from
stretcher to stretcher or from stretcher to
hospital bed.

It is very easy to apply the stretchers restraint
system due to the easy access slits at the
shoulder level on the Uni Easy Transfer
mattress.

For patient safety the mattress should be
secured on the stretcher with the stretchers
restraint system.

It is hygenic as it is fully welded and easy to
clean with most common disinfectants.